



## Infinite® 200 PRO: detection solutions that grow with your ideas

Immediate access to all wavelengths in an affordable, scalable detection family – with the patented NanoQuant Plate™ and Gas Control Module (GCM™)



NanoQuant Plate™  
compatible



## Infinite 200 PRO – New horizons in microplate detection











**Building on the success of the Infinite 200 series, Tecan has developed the Infinite 200 PRO, with new enhancements that cater to the needs of today's scientists. The Infinite 200 PRO offers flexible, scalable detection solutions for a wide range of assays, using monochromator- and filter-based technologies.**

### Access to a full range of leading detection methods

Infinite 200 PRO can provide a full range of leading detection methods in one easy-to-use modular instrument. Users can select from modules listed in the table below to create a perfect reader for their needs, without losing the option to upgrade in the future.





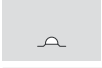




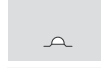

#### Infinite M200 PRO – Monochromator

The Quad4 Monochromators™ of the Infinite M200 PRO provides exceptional sensitivity, and allows the user to select any wavelength from UV to NIR, and to perform absorbance, excitation and emission scans. Users can access all wavelengths, and change from top to bottom reading, for easy measurement of multiplexed assays at the touch of a mouse click – no manual hardware changes are required.

	Fluorescence intensity top reading including TRF, with automated z-adjustment (optional background corrected z-focusing)
	Enhanced fluorescence intensity bottom reading with OR (Optimal Read) function, including TRF
	Spectrally enhanced photomultiplier tube
	Absorbance
	Photon counting luminescence, including dual color luminescence
	Cuvette port for absorbance
	Temperature control
	Injectors
	NanoQuant Plate
	Gas control module (GCM)

#### Infinite F200 PRO – Filter

The Infinite F200 PRO uses a patented intelligent filter slide system, with integrated lifetime monitoring of the filters, and its fluorescence polarization module is perfectly equipped for binding studies in homogenous mix and read assays. A dichroic filter allows TR-FRET applications, and the filter modules offer a cost-efficient solution for routine applications at fixed wavelengths.


	Fluorescence intensity top reading including TRF
	Enhanced fluorescence intensity bottom reading with OR (Optimal Read) function, including TRF
	TR-FRET/HTRF®
	Spectrally enhanced photomultiplier tube
	Absorbance
	Photon counting luminescence, including dual color luminescence
	Fluorescence polarization
	Temperature control
	Injectors
	NanoQuant Plate
	Gas control module (GCM)

## Select your new application, customize your detection device and perform your measurements quickly and easily

### Broadly applicable modular detection solutions to widen application capabilities

Detection is at the heart of biopharmaceutical and diagnostic assay measurements. In today's rapidly changing application environment the Infinite 200 PRO's modular, cost-effective design permits fast wavelength selection, and offers a host of easily upgradeable detection options.

The Infinite 200 PRO has been developed to deliver accuracy and performance in a format that allows you to build a versatile detection system to match your changing application needs. With the Quad4 Monochromators-based Infinite M200 PRO and filter-based Infinite F200 PRO detection options, the reader offers up to eight detection modes for sample measurements in 6- to 384-well plates, PCR plates or cuvettes. Three sets of advanced optics and three high performance detectors – optimized for the requirements of fluorescence, luminescence and absorbance reading – allow uncompromised performance in all detection modes.

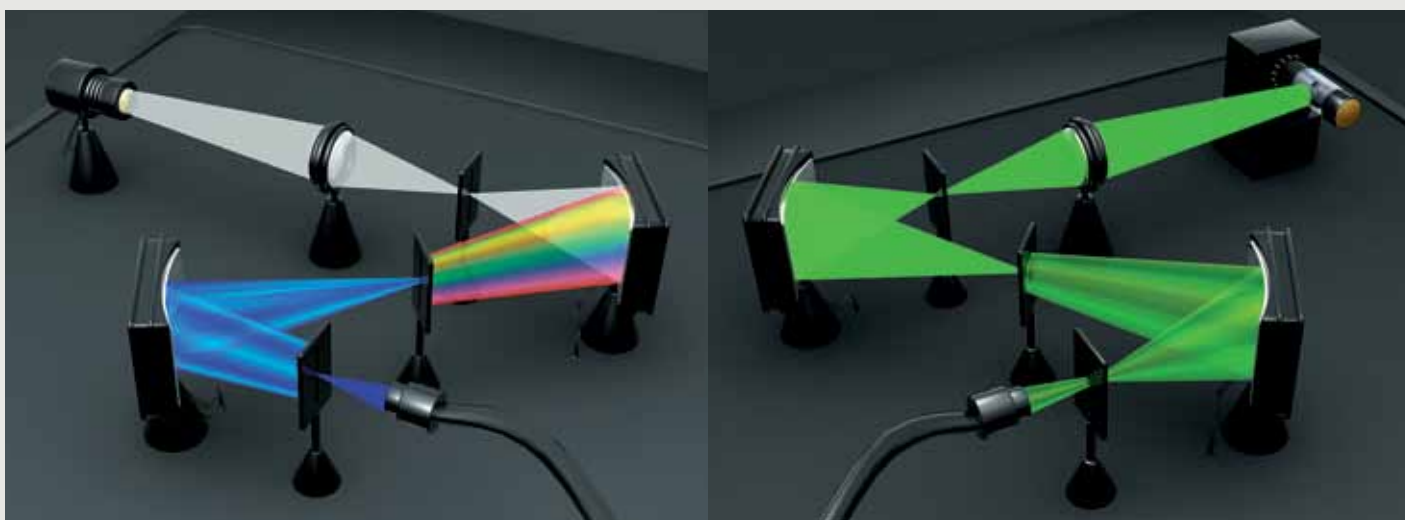

The Quad4 Monochromators technology makes use of a double monochromator on both the excitation and emission side. The picture above outlines the double monochromator system architecture on the excitation (left picture) and the emission (right picture).



Tecan's filter slide with patented system for monitoring of filter lifetime.

#### The Infinite 200 PRO offers unlimited flexibility for a wide range of biological assays and measurements including:

- DNA/RNA quantification
- protein quantification
- ion channel studies
- ion flux studies
- calcium ion detection
- reporter gene and gene expression assays
- cell viability and toxicity assays
- cell-based assays
- binding studies
- enzyme assays
- ELISA
- immunoassays
- fluorescence and luminescence applications
- TR-FRET/HTRF applications

## Various modules are available with the Infinite M200 PRO and Infinite F200 PRO Extended wavelength range and enhanced sensitivity

**A spectrally enhanced photomultiplier tube extends emission wavelength range from 330 – 600 nm to 280 – 850 nm, allowing the use of red-shifted dyes and minimizing interference caused by autofluorescence. A UV Si photodiode absorbance detector provides excellent sensitivity for the wavelength range of 230 – 1,000 nm, even at low concentrations.**

### Superior performance in absorbance for low sample volumes

The instrument's improved wavelength accuracy for 260/280 nm absorbance measurements allows high sensitivity determination of DNA or RNA concentration. Up to 16 samples with volumes as low as 2 µl can be measured simultaneously with Tecan's patented (EP2045015) NanoQuant Plate. This highly precise measurement tool uses a separate quartz optic for each sample, and requires no additional plate calibration.



### Gas Control Module

The patent pending Gas Control Module (GCM) for the Infinite 200 PRO offers a comprehensive solution for a variety of cell-based applications in this versatile multimode reader. Two independent gas inlets allow the automated control of CO<sub>2</sub> and O<sub>2</sub> concentration inside the reader chamber and help to maintain stable culture conditions during prolonged experiments and allow assays to be performed under anaerobic or physiological conditions. Maintaining the optimal CO<sub>2</sub> concentration within the measurement chamber helps stabilize pH and medium conditions, while the independent control of O<sub>2</sub> concentration (oxygen reduction is achieved using N<sub>2</sub>) provides hypoxic or simulated in vitro growth conditions. Combining this with precise temperature control and efficient shaking, the Infinite 200 PRO makes cellbased assays more biologically relevant. In addition, the elimination of data gaps (e.g. overnight or on weekends) minimizes the number of repeated assays and leads to more consistent and reliable data than can be achieved manually.

The GCM's acoustic and visible warning system offers excellent measurement reliability. This system is able to detect if the gas pressure or flow changes dramatically during an experiment or quickly recognize if the target concentration deviates or will not be reached.

Altitude influences the atmospheric partial pressure of CO<sub>2</sub>, affecting the measured value. The GCM's unique altitude correction function compensates for this, ensuring precise,


stable measurement and regulation of gas concentration inside the reader chamber.

Tecan's impressive GCM allows the optimization the gas mixture within the reader, providing the perfect solution for experiments with mammalian cells, hypoxia assays, cell viability studies, invasion assays, ischemia or reperfusion studies and many more.




### Comprehensive format flexibility

The Infinite M200 PRO offers outstanding format flexibility, and can perform both fixed wavelength and scanning spectrophotometric measurements, using standard 1 × 1 cm cuvettes or low volume microcuvettes in an upright position. In addition, it is compatible with all standard microplate formats, from 6- to 384-wells, including low volume plates and Tecan's unique NanoQuant Plate.

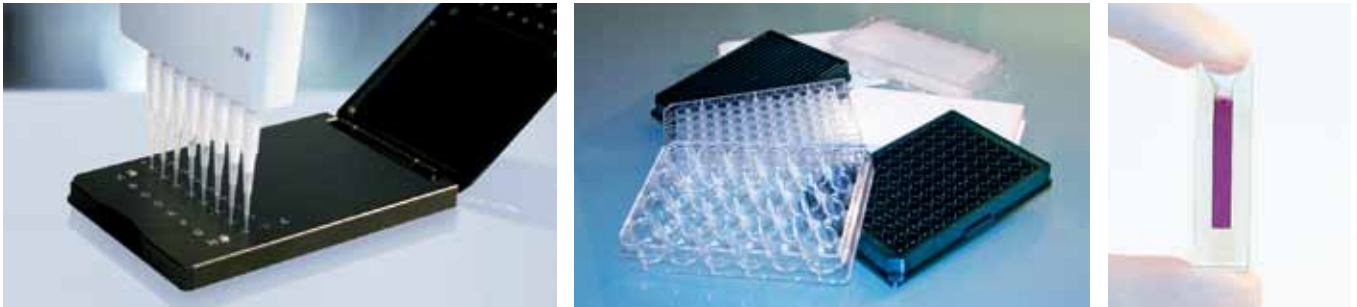


### Ready to go luminescence

The luminescence module is capable of reading dual-color luminescence assays, with a photon counting detector that can record even the lowest light levels from an assay, and an integrated set of luminescence filters enable BRET1 and BRET2 applications. The dynamic range for luminescence measurements has also been improved, helping the analysis of sets of samples with wide variation, without the need to adjust sample concentrations.

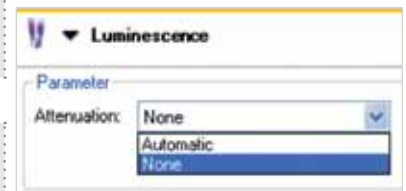


### Access to advanced assay systems

A dichroic mirror allows TR-FRET (HTRF) assays on the Infinite F200 PRO, and enhances detection limits for TRF Top Europium and FI Top Fluorescein measurements. This sophisticated system makes the Infinite F200 PRO an attractive and cost effective option for these demanding applications.

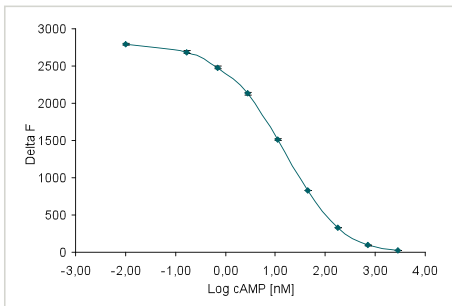


cAMP HiRange kit: The Delta F values obtained with the cAMP dilution series are inversely proportional to the cAMP concentration, resulting in the sigmoidal shape of the curve that is typical for competitive assays.

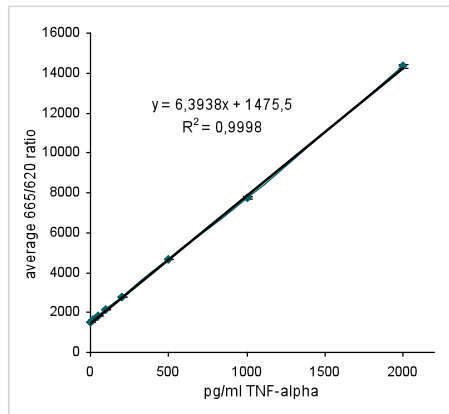

Human TNF-alpha kit: The measurement of a dilution series of the TNF-alpha standard shows a linear course ( $R^2 = 0.9998$ ) from 2000 to 20 pg/ml TNF-alpha.



### Automated, adjustable z-focus

Implementing assay miniaturization on the Infinite M200 PRO is helped by the automated, adjustable z-focus for FI Top measurements. Equally high sensitivity can be achieved for all plate formats, allowing the same high performance in low volume plates. This new feature, complete with background correction, is particularly suited to cell-based applications using autofluorescent growth media, providing automatic optimization of the signal-to-background ratio.

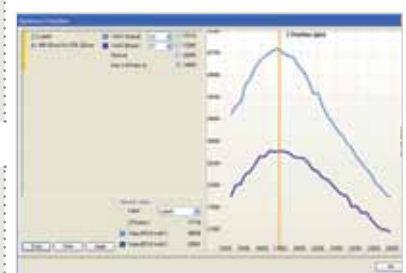

### Cell-based applications

The Infinite 200 PRO benefits from enhanced FI Bottom reading. Its special Optimal Read (OR) function has been designed specifically to optimize and improve cell-based measurements. Very low CVs, high intra- and inter-well reproducibility can be achieved when measuring adherent cells in microplates, offering increased sensitivity. The Infinite 200 PRO provides linear and orbital shaking – with adjustable amplitude in conjunction with frequency and duration – making it perfect for enzyme, bacterial and cell-based assays. The Infinite 200 PRO also allows temperature control for cellular and biochemical assays that require specific reaction temperatures, with top heating to avoid condensation in lidded plates, ensuring the best performance for covered MTP applications.

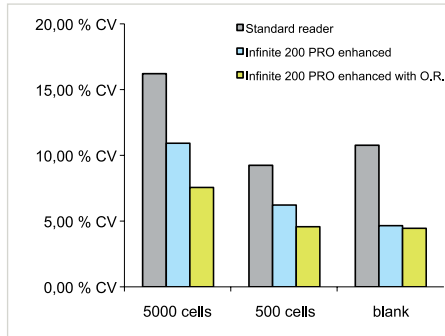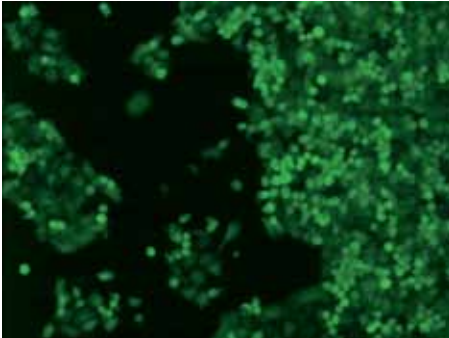

### Optimized injector module

The state-of-the-art injector module allows dispensing of up to two reagents per well, helping to replace a manual pipetting step or to trigger fast kinetic reactions in fluorescence, luminescence and absorbance modes. Its metal-free needles are ideal for ion studies, by preventing interference of metal ions in reactions. The injectors have variable volume and speed settings, and can be used together with the ratio mode to allow fast switching of wavelengths in a wide range of applications. The injector module has also been optimized for less wastage of substrates and buffers, with lower dead volumes for priming and the ability to tilt vessels, and its bulk reagent dispense function eliminates tedious pipetting steps for 6- to 384-well plates. Maintenance of the injectors is supported by easily accessible prime/wash buttons.



### MultiCheck™ – QC package for Infinite 200 PRO series

The Infinite 200 PRO has been designed to support users who need to meet Good Laboratory Practice standards. A MultiCheck QC plate, which includes installation and operational (IQ/OQ) checks and documentation, helps to ensure that all Infinite 200 PRO devices meet the standards needed for quality control laboratories, and satisfies the need to assure production standards in pharma and biotech settings.



### Built-in performance features


The new Plate In/Out button is another useful feature that has been introduced, in response to popular demand.




## Software designed for your workflow

Infinite 200 PRO users have complete access to intuitive software solutions that match their detection needs. The Infinite 200 PRO comes complete with i-control™ software interface that allows the user to define the workflow for each application.


Each workflow can be easily created by dragging and dropping the processing steps into the assay protocol sequence. The application workflow is then visible to the user, and can be saved for future use. Data sets are easily managed and exported to Windows® compatible formats like Excel®. The i-control software includes an application-oriented tab for rapid DNA/RNA quantification in the NanoQuant Plate, and identifies dye incorporation by measuring nucleic acid labeling efficiency. For more advanced data processing, Tecan's proven Magellan™ software provides features that perfectly match the flexibility of the Infinite 200 PRO. Magellan Tracker is designed to meet 21 CFR Part 11 requirements for electronic records and signatures, in compliance with FDA regulations.

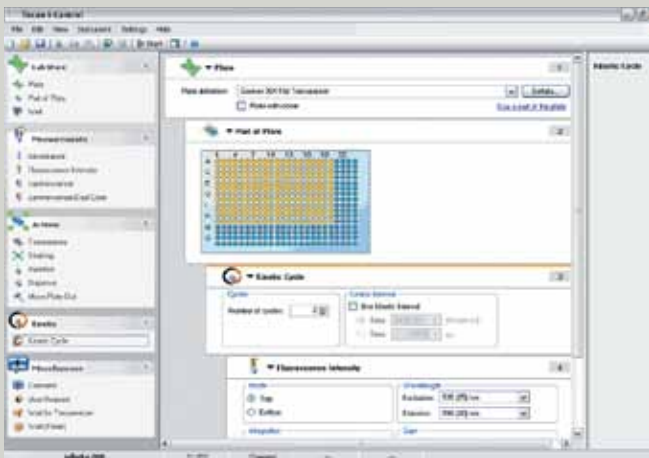


Workflow oriented i-control software supports complex assay protocols.

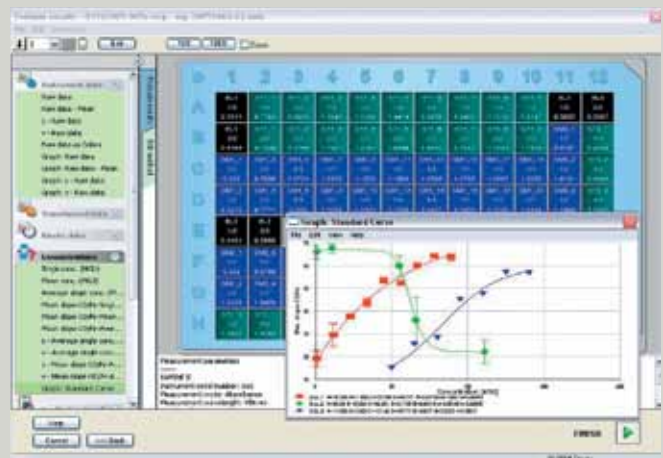

Magellan software allows easy presentation and evaluation of data from multiple experimental groups on a microplate.

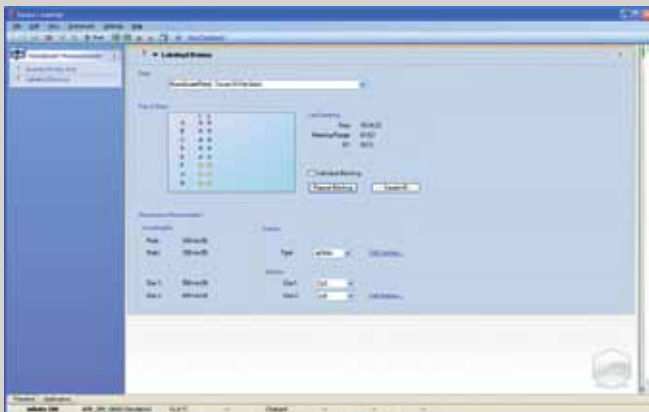

i-control application for nucleic acid quantification and measuring labeling efficiency.

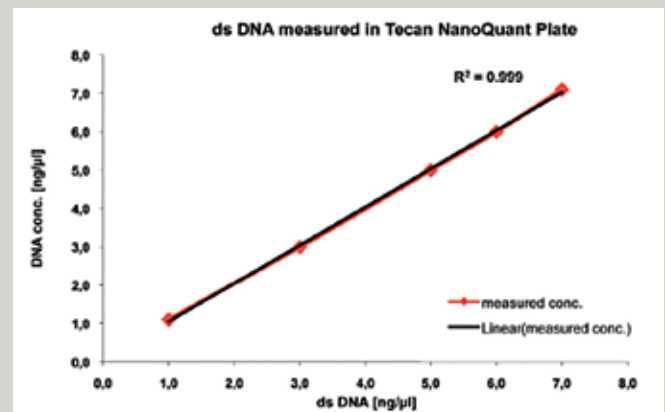

Linearity of low concentrated DNA dilutions.

### Highlights of new Magellan software in combination with the Infinite 200 PRO include:

- application-oriented workflow definition via drag-and-drop functionality
- wizard-guided application definition for intuitive operation, available in different languages
- easy conversion of data into results by Excel-style definition of transformations
- advanced spectra calculation package – the perfect partner for your Infinite M200 PRO reader
- convenient handling of dilution series and ICx calculations
- kinetic data analysis with calculation of slopes, onsets and enzyme kinetics
- pre-defined example files for a range of applications to help you get started immediately
- comprehensive plate library for fast selection of your favorite microplate.

# Infinite M200 PRO and F200 PRO – Typical performance values

<b>Light source</b>	UV Xenon flashlamp	
<b>Wavelength selection</b>	Quad4 Monochromators system (2 excitation and 2 emission monochromators)	
<b>Infinite M200 PRO</b>	Ex: < 5 nm for $\lambda \leq 315$ nm and < 9 nm for $\lambda > 315$ nm; Em: < 20 nm	
Bandwidth	<b>Absorbance</b>	<b>Fluorescence</b>
Wavelength accuracy	< $\pm 0,5$ nm for $\lambda > 315$ nm; < $\pm 0,3$ nm for $\lambda \leq 315$ nm	< $\pm 2$ nm for $\lambda > 315$ nm; < $\pm 1$ nm for $\lambda \leq 315$ nm
Wavelength reproducibility	< $\pm 0,5$ nm for $\lambda > 315$ nm; < $\pm 0,3$ nm for $\lambda \leq 315$ nm	< $\pm 1$ nm for $\lambda > 315$ nm; < $\pm 0,5$ nm for $\lambda \leq 315$ nm
<b>Infinite F200 PRO</b>	Up to 4 filter pairs per slide	
<b>Wavelength range</b>	Standard: Ex 230 – 600 nm, Em 330 – 600 nm	
Fluorescence intensity	Optional: Ex 230 – 850 nm, Em 280 – 850 nm	
Fluorescence polarization	Standard: Ex 300 – 600 nm; Em 330 – 600 nm	
Absorbance	Optional: Em 330 – 850 nm 230 – 1000 nm	
<b>Detectors</b>	Fluorescence – PMT, optional UV and red-sensitive	
	Absorbance – UV silicon photodiode	
	Luminescence – photon counting system with low dark current PMT	
<b>Plate formats</b>	6- to 384-well plates, cuvettes, NanoQuant Plate	
<b>Temperature control</b>	Ambient +5 °C up to 42 °C	
<b>Shaking</b>	Linear, orbital	
<b>Fluorescence sensitivity* values</b>	<b>Infinite F200 PRO</b>	<b>Infinite M200 PRO</b>
Fluorescence top reading*	85 amol / well (0,85 pM; 384-well plate)	170 amol / well (1,7 pM; 384-well plate)
Fluorescence bottom reading*	0,7 fmol / well (3,5 pM; 96-well plate)	1,2 fmol / well (6 pM; 96-well plate)
TRF**	2,8 amol / well (28 fM; 384-well plate)	90 amol / well (0,9 pM; 384-well plate)
FP*	< 4 mP standard deviation @ 1 nM Fluorescein	N / A
<b>Luminescence sensitivity values</b>	<b>Standard</b>	
Glow luminescence***	225 amol ATP / well (9 pM; low volume 384-well plate)	
Flash luminescence****	12 amol ATP / well (218 fM; 384-well plate)	
<b>Absorbance</b>		
Ratio accuracy 260 / 280 nm	$\pm 0,07$	
Precision @ 260 nm	< 0,2 %	
Accuracy @ 260 nm	< 0,5 %	
Measurement range	0 – 4 OD	
<b>Injectors</b>		
Pump speed	100 – 300 $\mu$ l/s	
Injection volume	selectable in 1 $\mu$ l increments; max. volume: 800 $\mu$ l per stroke	
Dead volume	100 $\mu$ l including pump back	
<b>Fastest Read Times</b>		
96 well plate	20 sec	
384 well plate	30 sec	
Wavelength Ex / Em-scan, 96 well plate	450 – 550 nm, 5 nm step	
	150 sec	

\* Detection limit for Fluorescein, \*\* Detection limit for Europium, \*\*\* Detection limit for ATP (144-041 ATP detection kit SL (BioThema)), \*\*\*\* Detection for ATP (ENLITE® Kit)

**Austria** +43 62 46 89 33 **Belgium** +32 15 42 13 19 **China** +86 21 2206 3206 **Denmark** +45 70 23 44 50 **France** +33 4 72 76 04 80 **Germany** +49 79 51 94 170  
**Italy** +39 02 92 44 790 **Japan** +81 44 556 73 11 **Netherlands** +31 18 34 48 174 **Singapore** +65 644 41 886 **Spain** +34 93 490 01 74 **Sweden** +46 31 75 44 000  
**Switzerland** +41 44 922 89 22 **UK** +44 118 9300 300 **USA** +1 919 361 5200 **Other countries** +41 44 922 8125

Tecan Group Ltd. makes every effort to include accurate and up-to-date information within this publication, however, it is possible that omissions or errors might have occurred. Tecan Group Ltd. cannot, therefore, make any representations or warranties, expressed or implied, as to the accuracy or completeness of the information provided in this publication. Changes in this publication can be made at any time without notice. All mentioned trademarks are protected by law. For technical details and detailed procedures of the specifications provided in this document please contact your Tecan representative. This brochure may contain reference to applications and products which are not available in all markets. Please check with your local sales representative.

Tecan and Infinite are registered trademarks and Magellan, NanoQuant Plate, Quad4 Monochromators, MultiCheck, GCM and i-control are trademarks of Tecan Group Ltd. Männedorf, Switzerland. Windows and Excel are registered trademarks of Microsoft Corporation in the United States and other countries. HTRF is a registered trademark of Cisbio bioassays, France. The fluorescence ratio associated with the HTRF readout is a correction method developed by Cisbio international and covered by the US patent 5,527,684 and its foreign equivalents, for which Cisbio international has granted a license to Tecan Group Ltd. Its application is strictly limited to the use of HTRF reagents and technology, excluding any other TR-FRET technologies. Transcreeper HTS Assay Platform is a patented technology of BellBrook Labs, LLC USA. Transcreeper is a registered trademark of BellBrook Labs, LLC USA. ENLITE Kit is a registered trademark and DLReady is a trademark of Promega Corporation, USA. LanthaScreen is a registered trademark of Invitrogen (life technologies) Corporation, Carlsbad, USA.

© 2012, Tecan Trading AG, Switzerland, all rights reserved. For disclaimer and trademarks please visit [www.tecan.com](http://www.tecan.com)

[www.tecan.com/Infinite200PRO](http://www.tecan.com/Infinite200PRO)

